



### SETUP

1. Place mattress on bedframe, ensuring umbilical cord is at foot end and 'feet' symbol is facing up.
2. Attach mattress straps to moving parts of bed frame. Fig 1
3. Position control unit on end of bed (or a flat surface). Fig 2
4. Attach umbilical cord to control unit. Fig 3
5. Attach mains lead. Fig 3
6. Switch on.

Fig 1



Fig 2

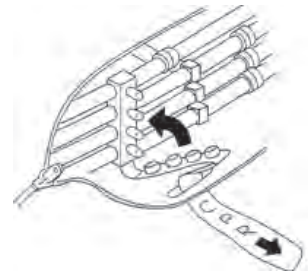


Fig 3










### CPR TAG


Ensure the CPR tag sealing connectors are pushed firmly onto the air pipes.




### OPERATION

- A  Power: Press for 2 secs to switch on/off.
- B  Alarm Mute: Mutes the audible alarm for 20mins.  
Refer to user guide for detailed instructions
- C  Alternating/Static: Switches between alternating and static. Reverts to alternating after 1hr
- D  Pressure Adjustment: Press '+/-' to increase or decrease the pressure setting.  
 (18mmHg - 60mmHg)

- E  Max Inflation: Rapidly inflates mattress to maximum static pressure. Reverts after 20mins.
- F  Lock/Unlock: Locks/unlocks the selected function. Mute remains operational. 'Max' and 'A/S' buttons still revert if selected.

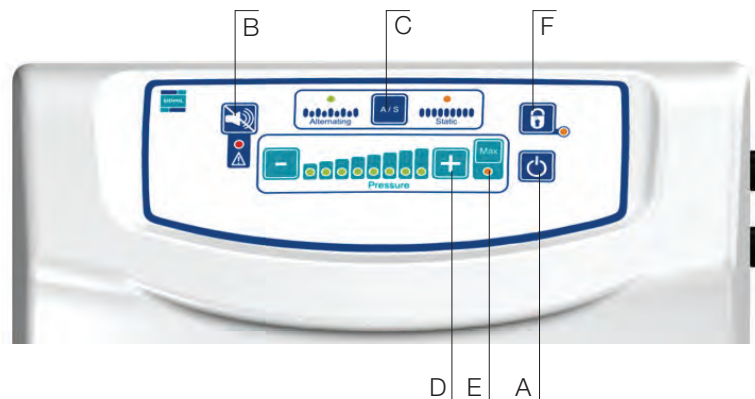

### Troubleshooting - Red alarm light will flash and audible alarm sound.

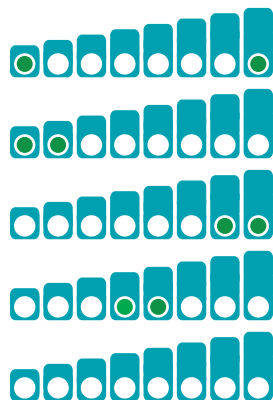

**Setup Failure:** Check CPR tag and umbilical cord are attached correctly. Check all air hoses/cells within mattress. Switch off/on.

**Pressure too Low:** Check CPR tag and umbilical cord are attached correctly. Check all air hoses/cells within mattress. Switch off/on.

**Pressure too High:** Disconnect umbilical cord to reduce pressure. Check for twists/kinks in air hoses. Switch off/on.

**Alternation Failure:** Disconnect umbilical cord to reduce pressure. Switch off/on.

**Power Failure:** Check mains cable is plugged in and switched on. Check fuse in plug and control unit.